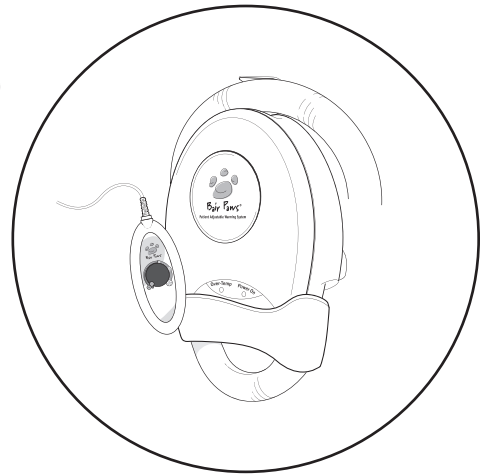
Bair Paws® Model 875 Warming Unit Operator's Manual

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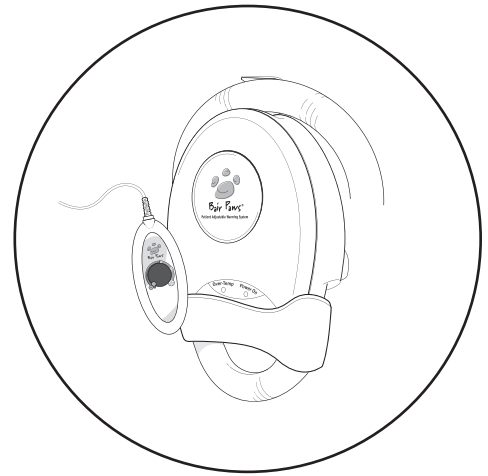
Table of Contents

Introduction to the Bair Paws® Patient Adjustable Warming System.....	1
Important Information about the Bair Paws Model 875 Warming Unit	2
Indications/Intended Use	2
Contraindication	2
Warnings	2
Cautions.....	3
Notices	3
Proper Use and Maintenance	3
Read Before Servicing Unit	3
Preparing the Bair Paws Model 875 Warming Unit for Use.....	4
Placing the Warming Unit on the Wall-mount Bracket	4
Mounting the Warming Unit on an IV Pole.....	5
Mounting the Warming Unit on a Bedrail	5
Mounting the Warming Unit with a Rail-mount to the Wall	6
Operating Information.....	7
Instructions for Use.....	7
Using the Hand-held Temperature Controller	8
Using the Holder for the Hand-held Temperature Controller	9
What to Do if the “Over-Temp” Indicator Light Illuminates and Audible Alarm Sounds	9
General Maintenance	10
Calibrating the Operating Temperatures	10
Replacing the Filter	10
Cleaning the Warming Unit.....	10
Technical Support and Customer Service	11
U.S. Customer Service.....	11
When You Call for Technical Support.....	11
Repair and Exchange.....	11
Technical Specifications	12
Definition of Symbols	14

Introduction to the Bair Paws® Patient Adjustable Warming System

The Bair Paws patient adjustable warming system consists of the Model 875 forced-air warming unit and disposable warming gowns. The Bair Paws Model 875 warming unit can provide comfort warming and prewarming to patients in preoperative settings. In postoperative settings the warming unit can provide comfort warming.

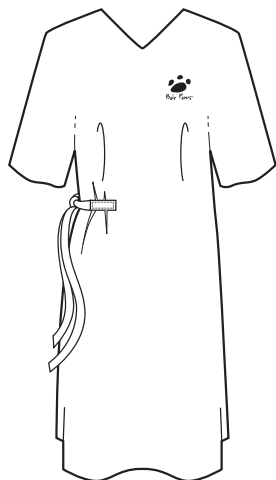
This manual includes operating instructions and unit specifications for the Model 875 warming unit. Please refer to the Instructions for Use for more information about the gowns.

Bair Paws Model 875 Warming Unit

The Bair Paws Model 875 warming unit has the following key elements: a blower, a heating element, and a hand-held temperature controller. The warming unit delivers warmed air through a hose that is connected to a port in a Bair Paws gown. The patient can adjust the air temperature and airflow using the temperature controller.

Bair Paws Warming Gowns

The Bair Paws single-use warming gown has an integral, channeled insert that delivers warm air through small perforations to warm the patient. Air channels and hose ports enable comfort and prewarming with a Bair Paws 800 series warming unit or clinical warming with a Bair Hugger® 500 or 700 series temperature management unit. The gown also has Velcro® fastener strips on each shoulder that provide easy access to the patient's arms and chest. The Bair Paws warming gowns are latex-free.

VELCRO is a registered trademark of Velcro Industries B.V.

Important Information about the Bair Paws Model 875 Warming Unit

Indications/Intended Use

The Model 875 warming unit is indicated for patient warming. This warming unit has been designed for use with the Bair Paws warming gowns in all preoperative and postoperative settings.

Contraindication

Do not apply heat to lower extremities during aortic cross-clamping. Thermal injury may occur if heat is applied to ischemic limbs.

Warnings

- Use only Bair Paws warming gowns with this warming unit. This warming unit has been designed to operate safely with Bair Paws warming gowns; use with other products may cause thermal injury. To the full extent permitted by law, the manufacturer and/or importer declines all responsibility for thermal injury resulting from the unit being used in conjunction with products other than Bair Paws warming gowns.
-  Do not warm patients with the warming unit's hose alone. Thermal injury may result. Always connect the hose to a Bair Paws warming gown before providing patient warming.
- Do not continue warming if the red Over-temp indicator light illuminates or an audible alarm sounds; thermal injury may result. Turn the warming unit off.
- Do not use a forced-air warming device over transdermal medications; increased drug delivery, patient death, or injury may occur.
- Do not allow the warming unit's hose to contact the patient's skin during patient warming. Thermal injury may result.
- Equipment not suitable for use in the presence of a flammable anaesthetic mixture with air or with oxygen or nitrous oxide.
- This warming unit is not intended for use in the operating room.
- Position the temperature controller cord and the hose away from the patient's neck or shoulders to avoid entanglement and/or injury.

Cautions

1. Do not initiate patient warming unless the warming unit is safely placed on a hard surface or securely mounted. Otherwise, injury may result.
2. Electrical shock hazard. Do not disassemble the warming unit unless you are a qualified service technician. There are electrically live parts within the unit when it is connected to a power source, even when the unit is in OFF or Standby mode.

Notices

1. The Bair Paws warming unit meets the international electronic interference requirements. If radio frequency interference with monitoring equipment occurs, connect the warming unit to a different power source.
2. Federal law (USA) restricts this device to sale by or on the order of a licensed healthcare professional.
3. To reliably ground this warming unit, only connect to receptacles marked “Hospital Only” or “Hospital Grade.”

Proper Use and Maintenance

Arizant Healthcare Inc. assumes no responsibility for the reliability, performance, or safety of the unit if:

- modifications or repairs are performed by unqualified personnel,
- the warming unit is used in a manner other than that described in the Operator’s Manual, or
- the warming unit is installed in an environment that does not meet the appropriate electrical requirements.

Read Before Servicing Unit

All repair, calibration, and servicing of the Model 875 warming unit must be performed by qualified, medical equipment service technicians who are familiar with good practice for medical device repair. If the warming unit does not require the manufacturer’s attention, Arizant Healthcare Inc. will ship replacement parts to your location. Perform all repairs and maintenance in accordance with the instructions provided with the replacement parts.

Preparing the Bair Paws Model 875 Warming Unit for Use

Before you begin using the warming unit, make sure that it is safely placed on a flat, hard surface such as a table, or securely mounted on a wall, IV pole, or bedrail.

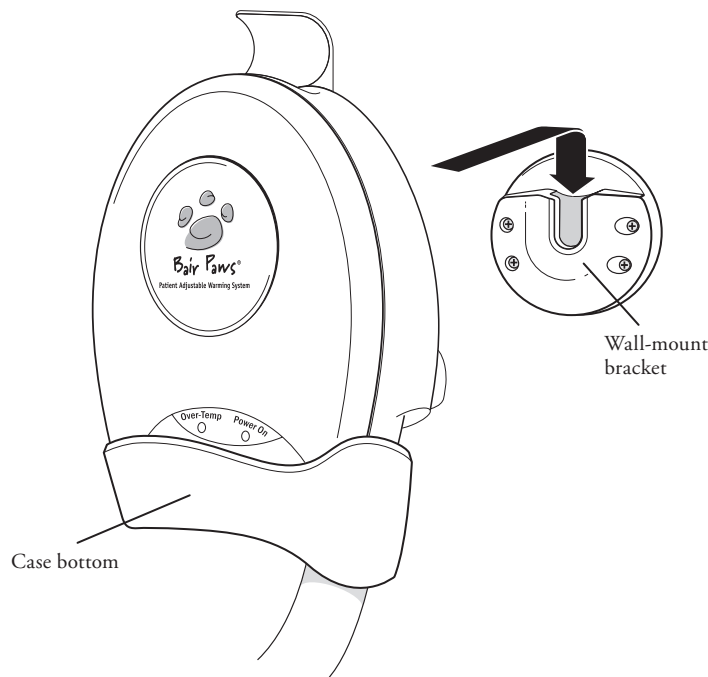
CAUTION: Do not place the warming unit on a soft or uneven surface, such as a bed; this may cause the air intake to become blocked, compromising the unit's performance.

A wall-mount bracket is shipped with every warming unit; the brackets for mounting the warming unit to an IV pole or a bedrail are available separately. Please contact your local sales representative or call Arizant Healthcare Customer Service at 1-800-733-7775 for more information about ordering mounting brackets.

The following instructions for mounting the warming unit assume that the brackets have already been attached to the wall or to the warming unit.

Placing the Warming Unit on the Wall-mount Bracket

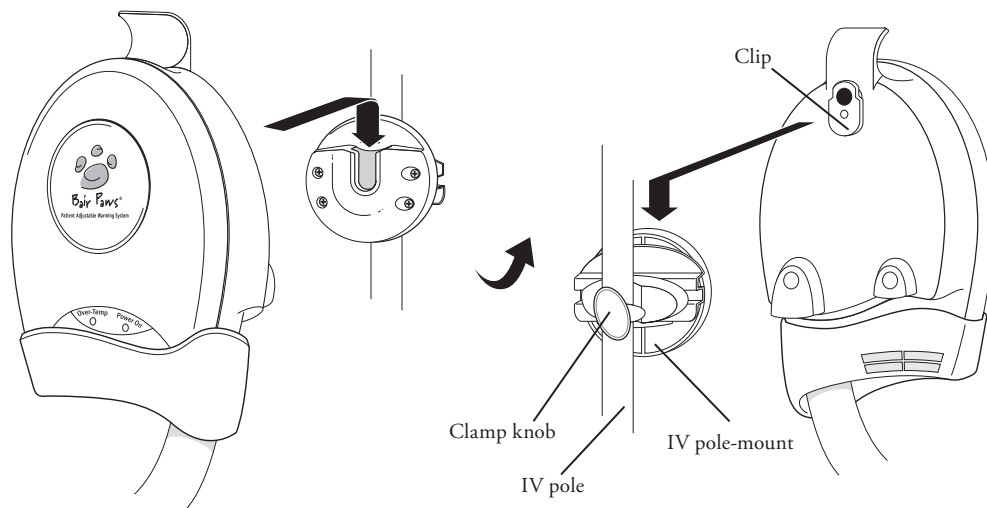
1. Snap the warming unit into position by sliding the clip on the back of the warming unit into the groove of the wall-mount bracket. The assembly will click when the clip is fully seated in the wall-mount bracket.
2. Lift the warming unit straight up to remove it from the wall-mount bracket.
NOTE: It may be necessary to gently strike the case bottom in an upward fashion with the palm of your hand to release the clip from the wall-mount bracket.

Mounting the Warming Unit on an IV Pole

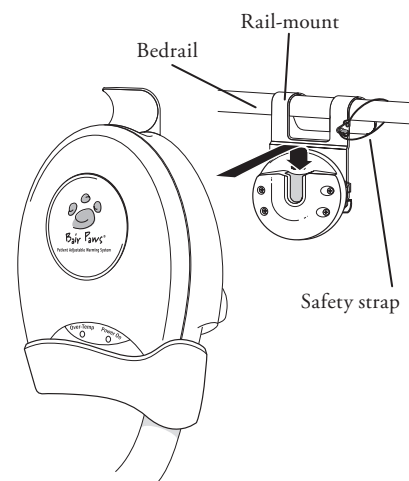
CAUTION: To prevent tipping, mount the Model 875 warming unit on an IV pole at a height that ensures stability. We recommend mounting the warming unit no higher than 44 in. (112 cm) from the floor on an IV pole with a minimum wheelbase radius of 14 in. (35.6 cm). Failure to do so may result in IV pole tipping, leading to potential catheter site trauma and other patient injuries.

1. Position the IV pole-mount at the desired height on the IV pole.
2. Turn the clamp knob on the IV pole-mount clockwise to tighten the clamp to the pole.
3. Snap the warming unit into position by sliding the clip on the back of the warming unit into the groove of the IV pole-mount. The assembly will click when the clip is fully seated in the IV pole-mount.

Mounting the Warming Unit on a Bedrail

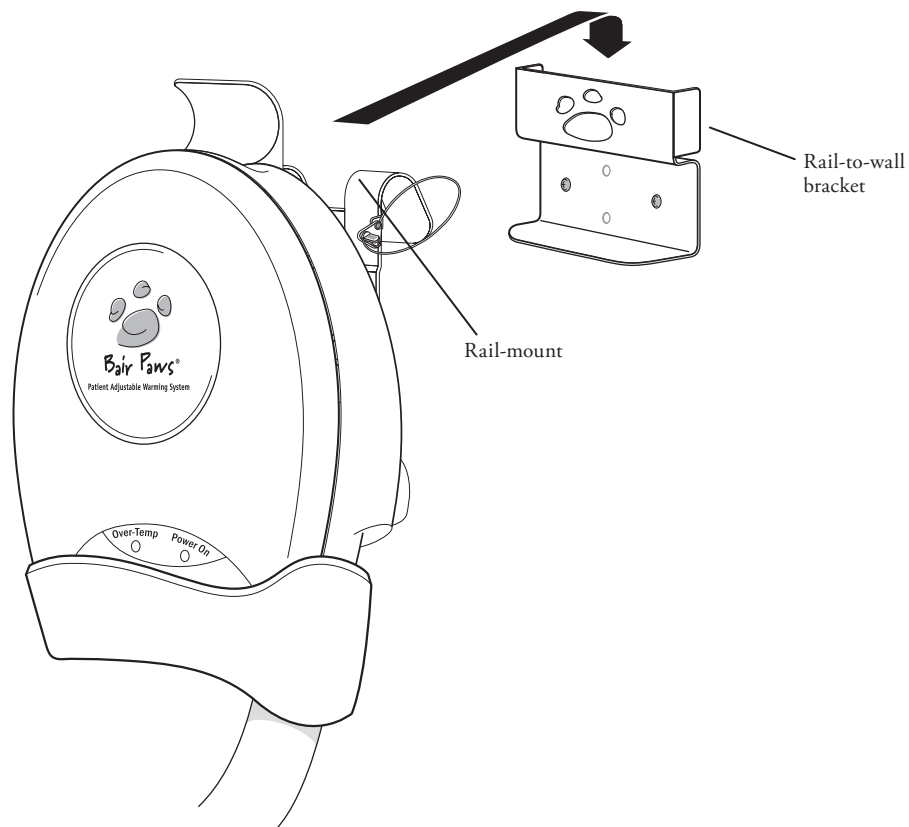
1. If the warming unit is not already connected to the rail-mount, slide the clip on the back of the warming unit into the groove of the rail-mount. The assembly will click when the clip is fully seated in the rail-mount.
2. Hang the Model 875 warming unit on a bedrail by the hooks on the rail-mount.
3. Wrap the safety strap around the bedrail and connect the end to the fastener on the strap. The safety strap will prevent the warming unit from falling if the unit is inadvertently dislodged from the bedrail.

Mounting the Warming Unit with a Rail-mount to the Wall

PRECAUTION: The rail-to-wall bracket is intended for use only with a Bair Paws 800 series warming unit.

1. If the warming unit is not already connected to the rail-mount, slide the clip on the back of the warming unit into the groove of the rail-mount. The assembly will click when the clip is fully seated in the rail-mount.
2. Hook the warming unit with the attached rail-mount to the rail-to-wall bracket.
3. Lift the warming unit straight up to remove it from the rail-to-wall bracket.

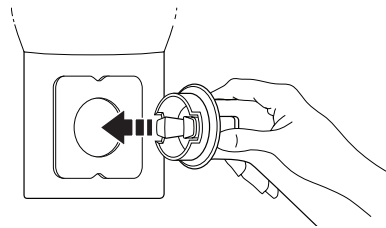
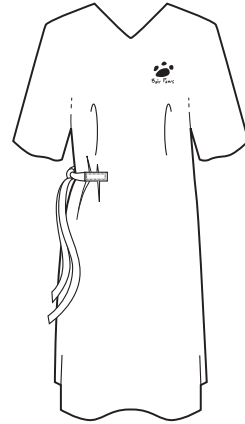
Operating Information

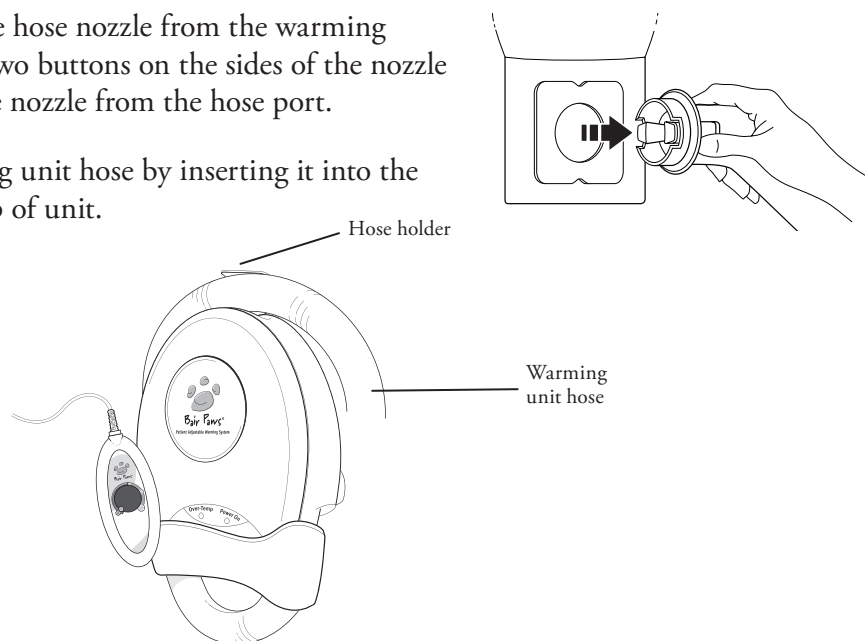
Instructions for Use

The instructions below describe how to use the Bair Paws Model 875 warming unit with a Bair Paws gown. For more information about the Bair Paws warming gowns, refer to the Instructions for Use that is shipped with the warming gowns.

1. Instruct your patient to put the warming gown on so the opening is in the back. The gown has ties on the side and back of the neck. The gown should be tied loosely to help the air circulate in the gown insert.
2. Connect the warming unit hose to the warming gown by performing the following steps.
 - a. Squeeze and hold the two buttons on the sides of the hose nozzle.
 - b. Insert the nozzle into the properly sized hose port in the warming gown. Refer to the Instructions for Use that is shipped with the warming gowns for hose port location.
 - c. Release the buttons to lock the hose nozzle in place.
3. Connect the warming unit to a proper power source.
4. Turn the knob on the temperature controller clockwise to turn the warming unit ON. When you turn the warming unit on, the heater and blower activate. Also, the “Power On” indicator light on the warming unit illuminates.
5. Allow the patient to adjust the temperature to a comfortable level. The temperature and airflow are adjusted by turning the knob on the temperature controller. The temperature adjustment varies from approximately ambient to $43\pm3^{\circ}\text{C}$.

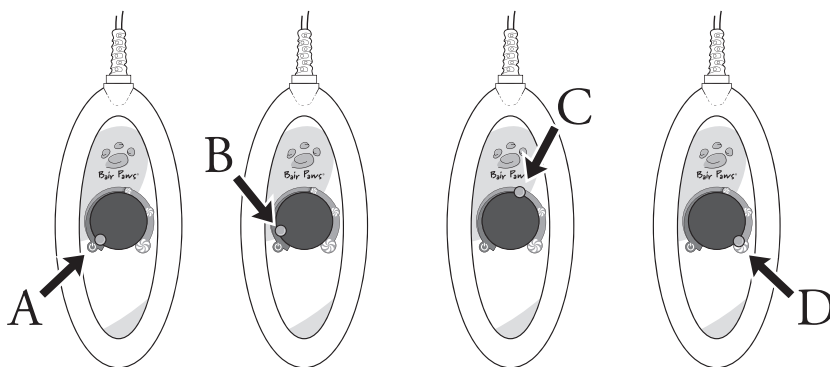
6. To disconnect the hose nozzle from the warming gown, press the two buttons on the sides of the nozzle and withdraw the nozzle from the hose port.
7. Store the warming unit hose by inserting it into the hose holder at top of unit.

Using the Hand-held Temperature Controller

After the hose has been connected to the Bair Paws gown, the patient can adjust the air temperature and airflow by turning the knob on the temperature controller. The following list explains how the unit will operate when the knob is placed in the specified position shown in the figure.

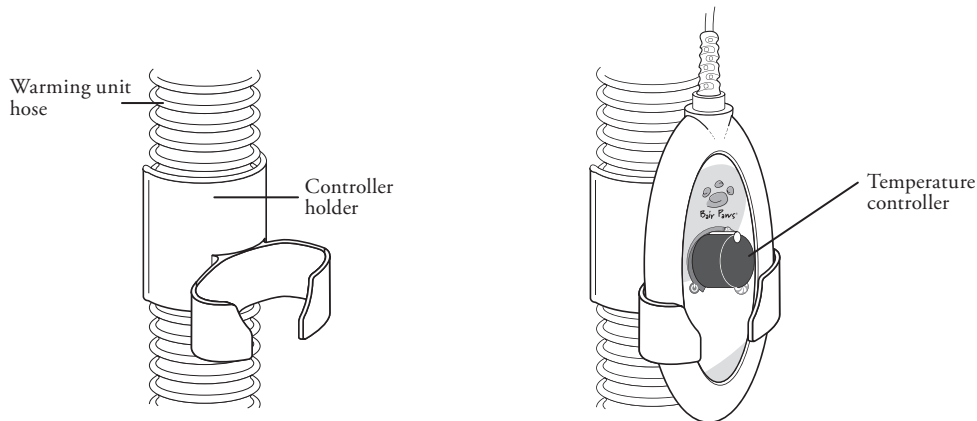
- When the knob is in position A, the warming unit is OFF.
- As the knob is rotated from position A to B, the knob clicks and the warming unit turns on LOW AIRFLOW and LOW HEAT.
- By turning the knob from position B to C, the air temperature increases from slightly above room temperature to approximately 41°C (106°F), while the blower remains on LOW AIRFLOW.
- As the knob rotates clockwise past position C, the blower's airflow increases.
- At position D, the warming unit is on HIGH HEAT and HIGH AIRFLOW. The air temperature will stabilize at 43±3°C (109±5.4°F).

Using the Holder for the Hand-held Temperature Controller

Store the hand-held temperature controller in the controller holder to help prevent cord clutter around the patient and accidental damage to the warming unit.

1. Press the holder onto the warming unit hose so it is within the patient's reach.
2. Slide the temperature controller into the holder. It should fit securely with the dial facing the patient.

What to Do if the “Over-Temp” Indicator Light Illuminates and Audible Alarm Sounds

If the “Over-Temp” indicator light illuminates, it means that the warming unit is delivering air at a higher temperature than the unit is calibrated to deliver. When an over-temperature condition occurs, the color indicator light illuminates, an audible alarm will sound, and the heater will automatically turn off (the blower will continue to run). You should perform the following steps:

1. Discontinue patient warming.
 - Disconnect the warming unit hose from the Bair Paws warming gown.
 - Turn the warming unit OFF.
2. Unplug the warming unit from the power source and wait 5 minutes.
3. Reconnect the warming unit to a proper power source.
4. Turn the warming unit ON and select the maximum airflow and temperature setting by turning the controller knob completely clockwise.
5. Allow the warming unit to run for at least 5 minutes.
 - If the “Over-Temp” light illuminates again, return the warming unit to Arizant Healthcare for service.
 - Otherwise, reconnect the warming unit hose to the Bair Paws gown to resume patient warming.

General Maintenance

Calibrating the Operating Temperatures

Calibrate the operating temperatures of the Model 875 warming unit (1) after performing a service procedure and (2) routinely every six months. Perform the calibration in a similar environment temperature to where the Model 875 will be used. You will need a Bair Paws Model 90055 temperature test kit to perform the calibration. Calibration instructions are included with the test kit.

Replacing the Filter

Replace the filter every 6 months or sooner if necessary. To order a replacement filter, contact Arizant Healthcare customer service using either of the telephone numbers listed on page 11. Instructions for replacing the filter will be included with the replacement filter.

Cleaning the Warming Unit

WARNING: Do not immerse any part of the warming unit in liquid while cleaning it. Moisture will damage the components, and thermal injury may result.

PRECAUTIONS:

- Do not use a dripping wet cloth to clean the warming unit. Moisture may seep into the electrical contacts and damage the components.
- Do not use harsh solvents to clean the warming unit.^a Solvents may damage the labels and other plastic parts.

METHOD

1. Disconnect the warming unit from the power source before cleaning.
2. Wipe the cabinet, the temperature controller, and the outside of the hose with a damp, soft cloth and a mild cleaning solution or antimicrobial spray.
3. Dry with a separate soft cloth.

^a WINDEX® multi-use liquid cleaner, ENVIROCID[®] disinfectant, and ALCONOX[®] industrial cleaner are three products that can be used to clean the equipment, as well as 3% household bleach, isopropanol, ethanol, ammonia or other phosphate-based glass cleaning solutions.

WINDEX is a registered trademark of S.C. Johnson & Son, Inc.
ENVIROCID[®] is a registered trademark of Metrex Research Corporation.
ALCONOX is a registered trademark of Alconox, Inc.

Technical Support and Customer Service

U.S. Customer Service

TEL:

800-733-7775

952-947-1200

FAX:

800-775-0002

952-947-1400

When You Call for Technical Support

Please be ready to give the Technical Support Representative the serial number of your Bair Paws warming unit. The serial number is located on the back of the unit.

Repair and Exchange

Call Arizant Healthcare customer service if your Model 875 warming unit requires service. A customer service representative will give you a Return Authorization (RA) number. Please use this RA number on all correspondence concerning your warming unit. Your customer service representative will also send a shipping carton to you at no charge upon request.

Technical Specifications

Physical Characteristics

DIMENSIONS OF WARMING UNIT	13 in. high x 4 in. deep x 7.7 in. wide 33 cm high x 10.2 cm deep x 19.6 cm wide
DIMENSIONS OF TEMPERATURE CONTROLLER	2.5 in. wide x 5.8 in. long 6.4 cm wide x 14.7 cm long
WEIGHT OF WARMING UNIT	7.0 lb; 3.2 kg.
MOUNTING OPTIONS	Wall mount bracket, IV pole clamp, and rail-mount bracket with safety strap.
HOSE	Detachable, flexible, and wipeable. 78 in. long x 1.5 in. wide; 198 cm long x 3.8 cm wide
FILTRATION SYSTEM	Dust filter included.
RECOMMENDED FILTER CHANGE	Change at least every 6 mo.

Temperature Characteristics

TEMPERATURE CONTROL	Electronically controlled using integrated circuit sensor.
HEAT GENERATED	1000 BTU/hr (average)
AVERAGE OPERATING TEMPERATURES AT THE END OF THE HOSE	User adjustable from: ambient to $43\pm3^{\circ}\text{C}$ ($109\pm5.4^{\circ}\text{F}$)
RECOMMENDED ENVIRONMENT TEMPERATURE RANGE	18° to 26°C (65° to 79°F)

Safety System Characteristics

THERMOSTAT	Independent electronic and heater (electromechanical)
OVERCURRENT PROTECTION	Dual input fused line.
SAFETY FEATURE	Over-temperature protection: color indicator light illuminates, audible alarm sounds, and heater shuts down (blower continues to run).
CERTIFICATIONS	IEC 60601-1: UL 60601-1; CAN/CSA-C22.2, No. 601.1; EN 60601-1-2; EN 55011; ASTM F2196-02

CLASSIFICATIONS

Classified under EN 60601-1 Guidelines (and other national versions of the Guidelines) as Class I, Type BF, Ordinary equipment, Continuous operation. Not suitable for use in the presence of flammable anesthetic mixtures with air or with oxygen or nitrous oxide. Classified by Underwriters Laboratories Inc. with respect to electric shock, fire, and mechanical hazards only, in accordance with UL 60601-1, ASTM F2196-02, and Canadian/CSA C22.2, No. 601.1. Classified under the Medical Device Directive (93/42/EEC) as a Class IIb device.

Electrical Characteristics

BLOWER MOTOR	Airflow: 7-13 cfm (3.3-6.1 $\frac{1}{s}$)
POWER CONSUMPTION	Peak: 600W Average: 340W
LEAKAGE CURRENT	Meets IEC 60601-1 and UL 60601-1 requirements
HEATING ELEMENT	480W Resistive
POWER CORD	15-foot, SJT, 3 cond., 15A 4.6 m, HAR, 3 cond., 10A
DEVICE RATINGS	110-120 VAC, 50/60 Hz, 4.6A 220-240 VAC, 50/60 Hz, 2.8 A
FUSES	6.0A (120VAC); 4.0A (220-230VAC)
TEMPERATURE CONTROLLER CORD	32" from hose collar, 4 cond., Max. voltage: 5V

Definition of Symbols

The following symbols may appear on the product, on the exterior packaging, or in the product labeling.

ON/STANDBY

Increase temperature and airflow

Temperature Control

Equipotentiality plug (ground)

Fuse

Attention, consult accompanying documents

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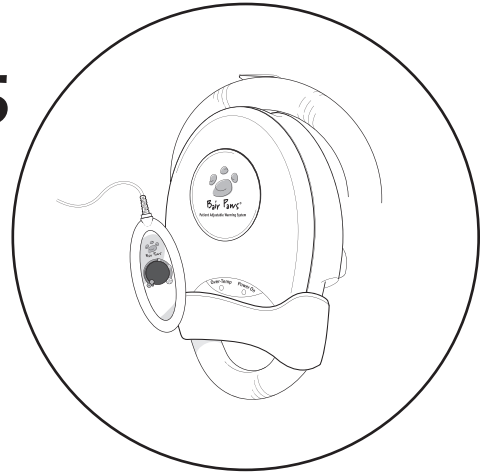
^a WINDEX® multi-use liquid cleaner, ENVIROCIDE® disinfectant, and ALCONOX® industrial cleaner are three products that can be used to clean the equipment, as well as 3% household bleach, isopropanol, ethanol, ammonia or other phosphate-based glass cleaning solutions.

WINDEX is a registered trademark of S.C. Johnson & Son, Inc.
ENVIROCIDE is a registered trademark of Metrex Research Corporation.
ALCONOX is a registered trademark of Alconox, Inc.

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^a WINDEX® multi-use liquid cleaner, ENVIROCID® disinfectant, and ALCONOX® industrial cleaner are three products that can be used to clean the equipment, as well as 3% household bleach, isopropanol, ethanol, ammonia or other phosphate-based glass cleaning solutions.

WINDEX is a registered trademark of S.C. Johnson & Son, Inc.
ENVIROCID is a registered trademark of Metrex Research Corporation.
ALCONOX is a registered trademark of Alconox, Inc.

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Authorized Representative in the European Community (as defined
in Article 14 of the Medical Device Directive: 93/42/EEC):
Arizant UK Ltd., Calder Island Way, Wakefield, WF2 7AW
TEL (44) 1924 200550 • FAX (44) 1924 200518

Bair Hugger, Bair Paws, Arizant, Arizant Healthcare, bright
ideas that work, and the Bair Paws and Arizant logos are
trademarks of Arizant Healthcare Inc., registered or pending in
the U.S. Patent & Trademark Office and in other countries.

May be covered by one or more of the following patents: U.S.
Patent 7,226,454; 6,876,884; D485,338. Other patents pending.

201337B 05/08

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©2008 Arizant Healthcare Inc. All rights reserved